



REVISTA ARGENTINA DE MICROBIOLOGÍA

www.elsevier.es/ram



ORIGINAL ARTICLE

Improvement of bovine pestiviral diagnosis by the development of a cost-effective method for detecting viral RNA in fresh specimens and samples spotted in filter papers

Paula M. Favaro^{a,*}, Ana Molineri^b, Maria J. Dus Santos^{c,d}, Luis F. Calvinho^{a,b},
Andrea Pecora^{c,d}

^a Facultad de Ciencias Veterinarias, Universidad Nacional del Litoral, Santa Fe, Argentina

^b Instituto de Investigación de la Cadena Láctea (INTA-CONICET), Estación Experimental Agropecuaria Rafaela, Santa Fe, Argentina

^c Instituto de Virología e Innovaciones Tecnológicas (IVIT) CICVyA, INTA Castelar, Buenos Aires, Argentina

^d Consejo Nacional de Investigaciones Científicas y Técnicas, CONICET, Buenos Aires, Argentina

Received 16 March 2022; accepted 18 October 2022

KEYWORDS

Bovine pestiviruses;
Diagnosis;
Filter papers;
RT-qPCR

Abstract Bovine pestiviruses are the causative agents of bovine viral diarrhea, a disease that causes severe economic losses in cattle. The aim of this study was to improve their diagnosis by developing a RT-qPCR to detect bovine pestiviruses A, B and H; and to set up a protocol for collecting, shipping and preserving bovine pestiviral RNA on filter papers. The developed RT-qPCR showed high sensitivity in detecting these viruses in different matrices: viral stocks, semen and serum samples. With regard to the possibility of using the technique to test serum pools, it was possible to identify a positive serum sample within a pool containing 30 sera. In addition to evaluating the qPCR from fresh samples, the use of filter papers to sow bovine samples was analyzed. The sampling method on two different filter papers using bovine blood drops was a useful alternative for diagnostic purposes and allowed to preserve pestiviral RNA for up to 12 months under refrigeration.

© 2022 Asociación Argentina de Microbiología. Published by Elsevier España, S.L.U. This is an open access article under the CC BY-NC-ND license (<http://creativecommons.org/licenses/by-nc-nd/4.0/>).

* Corresponding author.

E-mail address: pfavaro@fcv.unl.edu.ar (P.M. Favaro).

<https://doi.org/10.1016/j.ram.2022.10.002>

0325-7541/© 2022 Asociación Argentina de Microbiología. Published by Elsevier España, S.L.U. This is an open access article under the CC BY-NC-ND license (<http://creativecommons.org/licenses/by-nc-nd/4.0/>).

Please cite this article as: P.M. Favaro, A. Molineri, M.J. Dus Santos et al., Improvement of bovine pestiviral diagnosis by the development of a cost-effective method for detecting viral RNA in fresh specimens and samples spotted in filter papers, Revista Argentina de Microbiología, <https://doi.org/10.1016/j.ram.2022.10.002>

PALABRAS CLAVE

Pestivirus bovinos;
Diagnóstico;
Papeles de filtro;
RT-qPCR

Optimización del diagnóstico de *Pestivirus* bovinos mediante el desarrollo de una metodología económica para detectar ARN viral en muestras frescas y fijadas en papeles de filtro

Resumen Los *Pestivirus* bovinos son los agentes causales de la diarrea viral bovina, una enfermedad que genera importantes pérdidas económicas en el ganado vacuno. El objetivo de este trabajo fue mejorar su diagnóstico mediante el desarrollo de una RT-qPCR para detectar los *Pestivirus* bovinos A, B y H y diseñar un protocolo de recolección, envío y conservación de ARN viral en papeles de filtro. La RT-qPCR desarrollada demostró alta sensibilidad en la detección de estos virus en diferentes matrices: stock viral, suero y semen. Respecto de la posibilidad de usar la técnica para testear *pools* de suero, fue posible identificar un suero positivo dentro de un *pool* compuesto por 30 sueros. Además de evaluar la qPCR en muestras frescas, se analizó el uso de papeles de filtro para sembrar muestras de bovinos. La metodología de toma de muestras en dos tipos de papeles de filtro usando gotas de sangre fue una alternativa útil para el diagnóstico y permitió conservar ARN viral por hasta 12 meses a temperaturas de refrigeración. © 2022 Asociación Argentina de Microbiología. Publicado por Elsevier España, S.L.U. Este es un artículo Open Access bajo la licencia CC BY-NC-ND (<http://creativecommons.org/licenses/by-nc-nd/4.0/>).

Introduction

Bovine pestiviruses cause severe economic losses in cattle worldwide. Pestiviruses belong to the family *Flaviviridae*, and are single-stranded RNA viruses, classified into the genus *Pestivirus*, also including *Pestivirus C* (*Classical swine fever virus*) and *Pestivirus D* (*Border disease virus*)¹⁵. Two species of bovine pestiviruses (also known as *Bovine viral diarrhea virus*) are known: *Pestivirus A* and *Pestivirus B* (previously known as BVDV-1 and BVDV-2, respectively)²⁹. Two decades ago, an emerging bovine pestivirus was reported²⁸, which was recently classified as *Pestivirus H* (HoBi-like pestivirus)²⁹. Each genotype is further classified into cytopathic (CP) and non-cytopathic (NCP) biotypes, according to their ability to cause damage in cell cultures¹¹.

Bovine pestiviruses can be transmitted to susceptible animals through all the secretions of an infected animal. Persistently infected animals (PI) play a fundamental role in the transmission of the virus since they are the main source of infection through their fluids to their cohabitants. They acquire the viral infection through a NCP strain that crosses the placenta during pregnancy. From days 45 to 125 of pregnancy the immune system of the fetus is not developed; therefore, it is incapable to elicit a proper immune response. Infection in this period of time leads to the development of a PI and seronegative calf³⁵. Thus, PI animals spread the virus through their secretions and excretions lifelong.

Pestiviruses A, B and *H* are the causative agents of bovine viral diarrhea (BVD), which includes a wide range of syndromes, such as diarrhea, respiratory disease, reproductive failures, embryonic mortality, abortions, malformations and hemorrhagic syndromes, depending on the age of the host and the viral strain. In addition, in PI animals, bovine pestiviruses can cause a severe syndrome called mucosal disease when they are re-infected with a CP strain of the virus¹².

Identification and elimination of PI animals is essential for the control of BVD in infected herds. Furthermore,

application of biosecurity measures, such as testing all animals arriving at the herd, testing semen, avoiding the contact with neighboring herds, controlling other species that could host the virus, among others, is necessary to prevent the reintroduction of pestiviruses in the herd³⁶. In some European countries, eradication of bovine pestiviruses was accomplished by applying the mentioned biosecurity measures, including or not vaccination programs¹⁷.

Detection of PI animals requires the use of direct diagnostic techniques. Among the available assays, viral isolation is the gold standard test; however, it is expensive and takes several days to yield results, hence, being currently disused in most diagnostic laboratories. Other techniques, such as enzyme-linked immunosorbent assays, direct immunofluorescence and immunohistochemistry, can be used to detect viral antigens. As an alternative, reverse transcription-polymerase chain reaction (RT-PCR) and quantitative RT-PCR (RT-qPCR) detect viral RNA with high sensitivity^{31,32}. qPCR has several advantages, such as fastness, and the possibility to pool samples from different animals⁷. Once the pools have been checked by RT-qPCR, if one of them is positive, it has to be split, and each sample must be tested individually. Compared with individual testing, the chance to test several animals as a pool is a great advantage in herds with low BVD prevalence, taking into account the saving of reagents and time consumed by technicians. Several strategies to pool samples according to the prevalence of BVD in herds have been proposed¹⁹. Currently, commercial in vitro diagnostic RT-qPCR kits are available to assure the detection of a bovine pestivirus positive animal within a pool of 50 samples; however, they may be inaccessible in some low-income countries because of their price and import constraints.

One of the critical points for molecular diagnostics is correct sample preservation. As bovine pestiviruses are very labile RNA viruses¹⁵, a cold chain must be ensured to avoid erroneous results in the RT-qPCR assay; therefore, samples must be taken and sent chilled to the laboratory as soon as possible. This is an inconvenient in herds located in places

far from the diagnostic laboratories. Blood samples should be maintained at 4 °C while semen should be sent in liquid nitrogen, which represents high costs and shipping logistics. A useful alternative to improve collecting and shipping samples is filter paper (cotton-based, cellulose paper), such as the Flinders Technology Associates cards (FTA® cards), especially designed to maintain nucleic acids stable at different temperatures and to protect them from UV rays²⁵. Since microorganisms present in samples spotted onto filter papers are immediately inactivated, the FTA system is completely safe and the cards could be shipped by postal mail, at room temperature (RT). FTA cards are widely used in different areas, such as biobanking, forensic genetics and molecular epidemiology research, since they preserve nucleic acids for long time periods. They have been effective for the detection of RNA viruses (*Foot-and-mouth disease virus*, avian influenza virus, *Porcine reproductive and respiratory syndrome virus*, *Rotavirus* and some arboviruses)^{2,9,10,14,33}. This technology is easy to use and non-invasive⁵, as a small volume of sample is spotted on a filter paper and venipuncture is not needed.

The aims of this study were to develop and validate a Sybr Green-based RT-qPCR test to assess detection and quantification of all bovine pestiviruses in serum and semen samples and to set up a strategy for collecting, shipping and preserving viral RNA on filter papers for BVD diagnosis and biobanking as a cost-effective method.

Materials and methods

Virus, plasmid and clinical samples

Viral strains

A CP *Pestivirus B* (isolate vs253) strain cultivated in Madin Darby bovine kidney cells (MDBK) available in the Virology Institute and Innovative Technologies (IVIT) from the Research Center of Veterinary Sciences (CICVyA), INTA Castelar, was used as viral stock to standardize the RT-qPCR and to spike specimens in order to mimic clinical samples.

Semen and serum samples

To detect RNA from bovine sera, blood samples without anticoagulant from healthy animals (pestivirus-negative by an in-house validated RT-PCR) were collected and sent to the laboratory at 4 °C. Sera was obtained after a five minute-centrifugation at 1500 × g and then dispensed in sterile 1.5 ml tubes. The sera tubes were stored at 4 °C or at –70 °C in case the sample was not immediately processed.

To assemble pools, aliquots of the same volume of each bovine serum were mixed in a 1.5 ml tube. The volume of each serum to be used in a pool was calculated according to the size of the pool, with a final volume of 200 µl.

To detect viral RNA from semen, 250 µl straws of extended semen from healthy animals (pestivirus-negative by an in-house validated RT-PCR) were shipped in liquid nitrogen and stored at –70 °C. All samples were thawed at RT before processing.

Sample preparation

The amount of virus used to spike sera and semen straws was calculated based on expected PI animal viral titers

(10⁴ TCID₅₀/ml) as determined by several authors^{3,16}. A titrated *Pestivirus B* stock was used to spike the samples.

Plasmid construction

For the qPCR standardization, a plasmid containing the pestiviral 5'untranslated region (5'UTR) of the viral genome was constructed. Partial 5'UTR of a local *Pestivirus A* isolate (25366) (obtained by RT-PCR) was inserted into a pGemT-easy vector (Promega, USA) following the manufacturer's recommendations. Competent *Escherichia coli* DH5-α strains were used for transformation with the ligation product. Transformation was carried out by thermal shock, following molecular biology standard methods²⁶. Plasmid was purified using the Wizard Plus SV Minipreps kit (Promega, USA) following the manufacturer's instructions. Plasmid identity was confirmed using restriction endonuclease digestion. Finally, the purification product was quantified by Nanodrop 2000 (Thermo Scientific, USA) and the number of copies was calculated.

Clinical samples

Twenty-eight (28) clinical samples collected from naturally infected animals that were confirmed positive by viral isolation in cell culture and subsequently direct immunofluorescence, were used to ensure viral RNA detection by the assay. The specimens consisted of 4 intestines, 7 fetuses, 3 nasal swabs, one ocular swab, 2 esophagus, one lung, 2 spleens and 8 blood samples (provided and diagnosed by Laboratorio Azul and Specialized Veterinary Diagnostic Service at INTA Balcarce, Argentina).

qPCR standardization

qPCR standardization was performed according to the Minimum Information for Publication of Quantitative Real-Time PCR Experiments (MIQE) guidelines⁴.

RNA extraction

RNA extraction was performed using the High Pure PCR Template Preparation kit (Roche, Switzerland) following the manufacturer's instructions. Positive and negative controls were included, consisting of 200 µl of *Pestivirus B* (isolate vs253) and 200 µl of elution buffer, respectively. The obtained RNA was immediately stored at –70 °C until processing.

Retro-transcription and qPCR reaction

RT-qPCR reactions were carried out simultaneously using the iTaq Universal SYBR Green One-Step Kit (Bio-Rad, USA). A 10 µl final volume mix was prepared with 5 µl of iTaq Universal SYBR Green reaction mix, 0.125 µl of iScript reverse transcriptase, 0.25 µl of forward and reverse primer (Genbiotech, Argentina), 2.625 µl of nuclease free water (Biodynamics, Argentina) and 2 µl of RNA obtained in a previous step. The cycling program consisted of 10 min at 50 °C for reverse transcription reaction, 1 min at 95 °C for polymerase activation and DNA denaturation, followed by 40 cycles of 15 s at 95 °C and 1 min at 60 °C, finalizing with a 0.5 °C increment from 65 to 95 °C every 2 s, using StepOne equipment (Applied Biosystems, USA). In each assay, positive and negative controls were included.

Primer selection and concentration adjustment

Primer pairs 189–389 were used to develop the qPCR¹⁸. Primer sequences were 5'-AGTCGTCAR(A/G)TGGTTCGAC-3' and 5'-TCCATGTGCCATGTACA-3' for 189 (forward primer) and 389 (reverse primer), respectively and amplify a 201bp sequence of the 5'UTR of bovine pestivirus genome. The expected amplicon was verified using BLAST (<https://blast.ncbi.nlm.nih.gov/Blast.cgi>) and aligned with reference NADL strain (GenBank accession number M31182) between positions 190 and 390 of the complete sequence.

To establish the most efficient dilution of primers to set up the test, 10 folded dilutions of the plasmid were prepared and several concentrations of primers were assessed (390 nM, 420 nM, 450 nM, 470 nM, 525 nM, 535 nM and 540 nM). Primer efficiency was confirmed for each experiment using the formula: $E = -1 + 10^{(-1/\text{slope})}$ ²², whereby a slope close to -3.32 indicates 100% efficiency. This ensured that PCR products were amplified at an efficient rate and experiments were comparable among them.

Analytical sensitivity and linear dynamic range of detection

Limit of detection (LOD) was calculated for the different matrices: plasmid, viral stock and spiked serum and semen. For that purpose: 10-fold dilutions were performed, from 10^9 to 10^0 copies of the plasmid (3 replicates) and from 10^5 to 10^0 copies of the viral stock, sera and semen (12 replicates). In addition, 3 replicates of pools of bovine sera containing 10, 20, 30, 40 and 50 individual samples (one serum of each pool was spiked to mimic a PI sample) were assessed in order to determine the efficiency in pooling sera. The linear range of the qPCR was established by a standard curve with 10-fold dilutions of the plasmid from 2.5 to 2 500 000 000 copies.

Assay precision

Replication experiments were carried out to test positive controls (plasmid and viral stock) and to compare cycle threshold (CT) results in order to ensure precise detection of the target sequence. CT values of 25 repetitions of viral stock and 10 repetitions of 5'UTR plasmid (processed in triplicate) were registered. Mean, standard deviation, weighted average variation coefficient (CVpp) and maximum random error allowed (CVA%) were calculated.

Analytical specificity

Bovine respiratory syncytial virus (BRSV), *Bovine parainfluenza-3* (PI-3), *Bovine herpesvirus-5* (BHV-5), *Bovine herpesvirus-1* (BHV-1), *Blue tongue virus* (BTV), *Bovine Rotavirus type A* (RVA) and *Mycoplasma* spp. were used to evaluate the specificity of the test. In addition, several bovine pestivirus strains available at the IVIT, INTA Castelar, were used to check the ability of the qPCR to identify all possible viral variants that were detected in Argentina. Viral isolates and reference strains included: 13 strains of *Pestivirus A* (4 from sub-genotype a and 9 from sub-genotype b), 4 strains of *Pestivirus B* and 5 strains of *Pestivirus H*.

Inter-laboratory assays

DNA aliquots of plasmid and negative controls, as well as reagents and instructions to run the test were sent to three

different laboratories to perform the qPCR. Two of them were located in the IVIT Institute and the other one was from a private veterinary diagnostic laboratory. Involved laboratories have trained technicians and separate areas to prepare PCR mixes and sow samples. The assay was also carried out in equipment available in these laboratories: iQ5 (Bio-Rad, USA), ABI7500 (Applied Biosystems, USA) and MyGo (IT-IS Life Science, Ireland) thermal cyclers.

Data analysis

Data analysis was carried out using StepOne Software version 2.3.

Filter paper sample processing

Storage of filter papers spotted with spiked blood

Two types of filter papers were used throughout the experiments: FTA cards (Whatman, UK) and cellulose 3mm chromatography paper (GE Healthcare, USA). Since FTA cards are designed for nucleic acid conservation over long time periods, they were used for conservation purposes. Conversely, chromatography paper, commonly used for other purposes and less expensive than the former, was evaluated for rapid laboratory tests in cases of suspected BVD in the field. This means, preserving RNA for shorter periods of time, allowing the collection and shipping of samples to the laboratory at room temperature.

Whole-blood samples were collected from healthy bovines, and then spiked with a titrated *Pestivirus B* viral stock as previously described so that there were 50 000 copies of RNA in a 50 μ l drop, according to PI animal viral titers.

Five repetitions of 50 μ l of spiked blood drops were spotted on chromatography paper using micropipettes and sterile tips. Then, the drops were dried at RT for 4h and finally stored at RT (15–25 °C), 4 °C, -20 °C and -70 °C for one and 7 days and 6, 12 and 18 months. The same methodology was performed on FTA cards: 5 drops of spiked blood were spotted for each temperature and stored for 12 and 18 months. The spiked filter papers were stored at different temperatures in zipped bags with desiccant to prevent humidity damage.

Once storage time expired, RNA was eluted. Briefly, four punches in each blood drop were cut with a 6mm puncher and placed in 1.5ml tubes containing 220 μ l of minimum essential medium (MEM). The tubes were then kept in a shaker for 30 min at RT and aliquots of 200 μ l were taken to proceed with RNA extraction and RT-qPCR described above including viral stock used as a positive control.

In order to prevent cross-contamination among samples and replicates, a protocol recommended in another study was followed. Briefly, the puncher was decontaminated between samples by dipping it in isopropyl alcohol for a minute. Later, the puncher was dried with absorbent paper and five punches were cut in a clean, new chromatography paper²⁷.

Limit of detection in filter papers

To assess the LOD in blood spotted on filter papers, samples were spiked with viral stock at several concentrations and spotted on filter papers. For this purpose, two

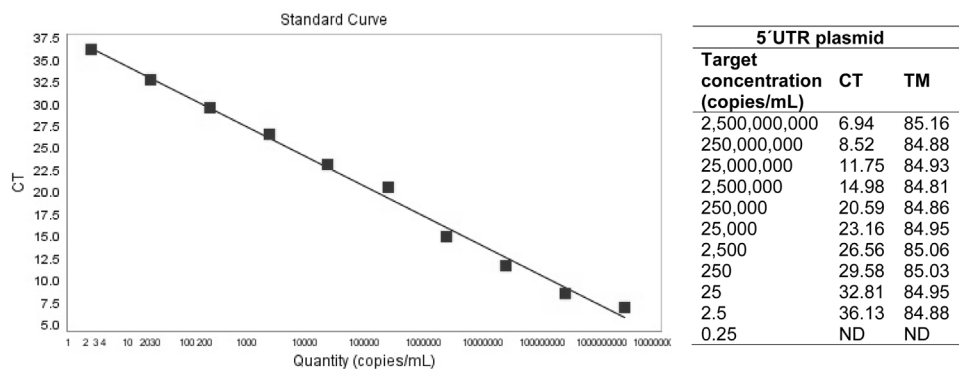


Figure 1 Left: Plasmid linear range of detection by the developed qPCR. Right: CT (cycle threshold) and TM (melting temperature in °C) values for 10-folded dilutions of the 5'UTR plasmid. ND: not detected.

replicates of 50 µl drops of blood containing 100, 1000, 5000, 10 000, 20 000, 40 000 and 80 000 viral particles were spotted on chromatography paper and later processed as mentioned above.

Recovering RNA from clinical samples

Blood drops of two PI animals were spotted on chromatography paper, dried up at RT for 4 h and then stored at 4 °C for 1 day and 6, 12 and 18 months. Once storage time expired, samples were processed as previously mentioned.

Statistical analysis

Aiming to compare the efficiency in recovering RNA between both types of papers at 12-month storage, a statistic analysis was conducted using the T-student and Fisher's tests. CT values and positive/negative results for each repetition stored at 4 °C and -20 °C were used for the T-Student and Fisher's tests, respectively.

Internal control to monitor nucleic acid extraction and PCR inhibitors

The efficiency of nucleic acid extraction and the presence of inhibitors for PCR amplification in semen and blood spotted on filter papers, fresh semen and tissue samples was checked by a qPCR assay using a set of primers that specifically amplify a fragment of the constitutive bovine GAPDH gene. For this purpose, GAPDH forward and reverse primer pairs²⁴ were used.

Results

RT-qPCR standardization

With regard to primer concentration adjustment, the best efficiency was observed in 470 nM (0.47 µM). The linear range of the PCR was from 2.5×10^9 to 2.5×10^1 DNA copies/ml, with an efficiency of 97.3%, a slope of -3.389 and regression coefficient of 0.994. Based on the calculations made on the plasmid, the test detected 2.5 copies (LOD) (Fig. 1). Regarding the qPCR sensitivity for viral stocks, the qPCR amplified up to five copies of viral RNA. Among the sera and semen matrices, the test was able to detect up

to 50 RNA copies (Table 1a). Furthermore, the RT-qPCR test was able to detect a positive serum in a pool containing 30 bovine sera (Table 1b).

In regard to specificity, the qPCR did not amplify unrelated microorganisms (BRSV, PI-3, BHV-1, BHV-5, BTv, RVA and *Mycoplasma* spp.). Furthermore, the test was able to detect all bovine pestivirus strains evaluated. Melting temperature (TM) values for each bovine pestivirus strain assayed is shown in Table 2, the specific range of TM went from 84.69 to 85.78 °C. The precision analysis of the 5'UTR plasmid and viral stocks showed a relatively low standard deviation and poor dispersion of CT results with respect to the mean (Table 3). Concerning the inter-laboratory assays, it was possible to detect the plasmid in the three different laboratories using several cyclers, with CT values ranging from 13.92 to 22.

All clinical specimens collected from naturally infected animals which were positive to conventional gold standard methods (viral isolation followed by direct immunofluorescence) were positive when analyzed through the developed RT-qPCR, with CT values ranging from 15 to 35 (Supplementary material 1).

Finally, the GAPDH gene, which was used to monitor nucleic acid extraction and possible presence of PCR inhibitors in blood and semen spotted on filter papers, fresh semen and tissue samples, were positive in all runs, with TM results ranging from 80 to 81 °C and homogeneous CT values (20 ± 3) for all specimens.

Pestivirus detection on filter papers

Filter papers spotted with spiked blood

As expected, when the storage period increased, the RNA recovering percentage among repetitions decreased. For instance, one and 7-day storage of blood drops on chromatography papers yielded a higher recovery percentage with respect to positive control CT when compared with 6 and 12-month storage. Regarding storage temperatures, in general, 4 °C, -20 °C and -70 °C were more effective in preserving RNA that RT, with a more marked trend as storage time increased (Table 4a).

When analyzing the samples stored for 12 months, results were similar in both types of papers (chromatography paper and FTA cards) at 4 °C and -20 °C storage temperatures

Table 1 (a) Number of copies detected by the qPCR in different matrices: virus stock (*Pestivirus B vs253*), serum and semen straws. (b) CT (cycle threshold) mean for three replicates of the different pool sizes evaluated.

| (a) Target concentration (copies/ml) | Viral stock | Serum | Semen straws |
|--------------------------------------|-------------|-------|--------------|
| 500 000 | + | + | + |
| 50 000 | + | + | + |
| 5000 | + | + | + |
| 500 | + | + | + |
| 50 | + | + | + |
| 5 | + | ND | ND |

| (b) Pool size (number of individual sera within the pool) | CT |
|---|------|
| 10 | 30.3 |
| 20 | 30.7 |
| 30 | 32.2 |
| 40 | ND |
| 50 | ND |

ND: not detected.

Results are expressed as 95% agreement between 12 replicates.

Table 2 Melting temperature (TM) of bovine pestiviral strains assayed by the RT-qPCR.

| Strain | TM (°C) |
|---|---------|
| Bovine Pestivirus A ^a Singer | 85.41 |
| Bovine Pestivirus A ^a MF120553 | 85.50 |
| Bovine Pestivirus A ^a 88625 | 85.78 |
| Bovine Pestivirus A ^a 71267 | 85.2 |
| Bovine Pestivirus A ^b 63588 | 85.08 |
| Bovine Pestivirus A ^b 90611 | 85.59 |
| Bovine Pestivirus A ^b MF120552 | 85.10 |
| Bovine Pestivirus A ^b 59473 | 85.32 |
| Bovine Pestivirus A ^b 83532 | 85.62 |
| Bovine Pestivirus A ^b 73611 | 84.81 |
| Bovine Pestivirus A ^b MF120594 | 85.22 |
| Bovine Pestivirus A ^b 78256 | 85.42 |
| Bovine Pestivirus A ^b 25366 | 85.32 |
| Bovine Pestivirus B 76/08 | 85.36 |
| Bovine Pestivirus B vs253 | 85.59 |
| Bovine Pestivirus B MF120585 | 85.44 |
| Bovine Pestivirus B MF120586 | 85.61 |
| Bovine Pestivirus H LV210185K/13 | 84.69 |
| Bovine Pestivirus H MK017821 | 84.80 |
| Bovine Pestivirus H MH992643 | 85.58 |
| Bovine Pestivirus H #48 | 85.1 |
| Bovine Pestivirus H #51 | 84.92 |

^a Pestivirus A (formerly called BVDV-1) sub-genotype a.

^b Pestivirus A (formerly called BVDV-1) sub-genotype b.

(Table 4b). No significant differences were identified between both types of papers at 12-month storage at 4 °C ($p=0.865$) and at -20 °C ($p=0.776$) using the T-Student test, neither at 4 °C ($p=0.500$) nor at -20 °C ($p=0.500$) using the Fisher's test.

Unfortunately, it was not possible to process the samples stored at -70 °C. At 18 months, results were negative for all storage temperature variables in both types of papers.

Limit of detection in chromatography paper

With regard to the LOD in blood samples, the qPCR was able to detect up to 1000 copies of viral RNA (Table 5).

Recovering RNA in papers from clinical samples

Blood drops from two PI animals spotted on chromatography paper and stored for one day and 6 and 12 months at 4 °C yielded positive results by the developed RT-qPCR, suggesting that after one year of storage at 4 °C, there was no significant degradation of the viral RNA on the chromatography paper for both PI animals. CT values for the samples for one day storage were 28.39 and 32.9; for 6-month storage were 23.93 and 33.75; finally, for 12-month storage were 29.93 and 27.56, respectively. On the contrary, results were negative at 18-month storage time.

Discussion

The developed RT-qPCR was able to detect the viral RNA in all bovine pestiviral strains tested and showed high specificity. A major advantage is the ability to detect *Pestivirus H* strains, which can be useful in areas where *Pestivirus H* are co-circulating with *Pestivirus A* and *B*. Furthermore, complete agreement was observed between RT-qPCR and viral isolation results of all clinical samples tested. In pooled sera, the test was able to detect a positive animal among 30 sera, which is a significant advantage when carrying out BVD control in low-prevalence herds since it allows to reduce the number of assays. However, in herds where BVDV prevalence is high, there is an increased chance that any single pool will test positive if the pool size is too large, requiring additional testing to identify the viremic individuals¹⁹. Thus, prevalence estimation should be considered when laboratories use the pooling-strategy, especially in regions with high BVD prevalence³⁰. Moreover, when the pool contains a large number of individual samples, the dilution effect can cause false negatives. This research has provided an in-house validated RT-qPCR that can replace commercially available kits, which represent a great trading cost for some countries,

Table 3 Mean, standard deviation, repetitions number (n), weighted average variation coefficient (CVpp) and maximum random error allowed (CVa%) of positive controls repetitions.

| | Mean | Standard deviation | n | CVpp | CVa% |
|--|-------|--------------------|----|-------|------|
| 5'UTR plasmid (2.5×10^3 copies) | 21.36 | 1.49 | 10 | 14.94 | 8.96 |
| <i>Pestivirus B</i> vs253 (5.62×10^6 copies) | 13.78 | 1.94 | 25 | 10.63 | 6.38 |

Table 4 Replicates cycle threshold (CT) means, number of positive determinations and recovery percentage at different time-points and temperatures for spiked blood spotted on filter papers. (a) From chromatography paper. (b) From FTA cards.

| | RT | | | 4°C | | | -20°C | | | -70°C | | |
|------------|---------|----|------|---------|----|------|---------|----|------|---------|----|------|
| | CT mean | PD | SD | CT mean | PD | SD | CT mean | PD | SD | CT mean | PD | SD |
| <i>(a)</i> | | | | | | | | | | | | |
| 1 day | 27.14 | 5 | 2.09 | 25.58 | 4 | 0.57 | 25.64 | 5 | 1.98 | 27.23 | 4 | 1.78 |
| Recovery % | 84.16 | | | 89.29 | | | 89.08 | | | 83.88 | | |
| 7 days | 28.77 | 5 | 4.67 | 28.67 | 5 | 0.99 | 28.47 | 5 | 1.12 | 26.26 | 5 | 2.12 |
| Recovery % | 79.39 | | | 79.67 | | | 80.22 | | | 86.98 | | |
| 6 months | 32.1 | 3 | 1.18 | 25.12 | 3 | 1.73 | 24.83 | 5 | 0.84 | 22.02 | 5 | 0.99 |
| Recovery % | 48.35 | | | 61.78 | | | 62.51 | | | 57.67 | | |
| 12 months | ND | 0 | - | 27 | 4 | 0.6 | 26.87 | 4 | 0.96 | | | NP |
| Recovery % | 0 | | | 67.81 | | | 68.14 | | | | | |
| <i>(b)</i> | | | | | | | | | | | | |
| 12 months | ND | 0 | - | 27.11 | 3 | 0.55 | 28.52 | 3 | 8.78 | | | NP |
| Recovery % | 0 | | | 67.81 | | | 64.20 | | | | | |

ND: not detected; NP: not processed.

Table 5 Values for spiked blood spotted on chromatography paper at several target concentrations.

| Target concentration (copies/ml) | Blood |
|----------------------------------|-------|
| 80 000 | + |
| 40 000 | + |
| 20 000 | + |
| 10 000 | + |
| 5 000 | + |
| 1 000 | + |
| 100 | ND |

ND: not detected.

Results are expressed as 95% agreement between 12 replicates.

and can be useful for veterinary diagnostic laboratories equipped with the infrastructure to perform qPCR assays.

Because of the SYBR Green-based detection system, the developed RT-qPCR presents some advantages: it is sensitive, easy to use and inexpensive. However, it is important to recognize the presence of any double-stranded DNA (e.g., primer-dimers) that can cause false positive results. This inconvenient can be sorted including an analysis after the amplification run, using the melting temperature of the expected amplicon, which allows the discrimination of the target amplicon to undesired products that can interfere with the results²³. In addition, using One-Step for retro-transcription and PCR reaction simultaneously reduces the

costs, time consumed and possibilities of contamination of the assay.

Diagnostic laboratories routinely use blood or serum for BVD analyses. In terms of practicability in the use of filter papers at the point of care, whole blood is more suitable than serum since it avoids the centrifugation step. Although semen is not used as a sample of choice to diagnose bovine pestiviruses, it is useful for the control of BVD in breeding herds^{20,21}; however, our group tested the same methodology using semen straw samples spotted on filter papers without achieving satisfactory results (Supplementary material 2).

Concerning the collection and conservation of bovine pestiviral RNA on filter paper, the literature is scarce. One article describes the use of several filter papers to preserve RNA up to six months in blood and sera samples from PI animals at RT, 4°C and -18°C³⁴. In this study, the authors evaluated four different types of papers: classical filter papers, Whatman paper No. 1, nitrocellulose membrane and HYBOND™-M nylon membrane. Results were similar for all papers tested except for classical filter papers, which yielded lower PCR products. However, the authors used an end-point RT-PCR assay and RNA isolation was carried out directly from the filter papers, avoiding the elution step. Unfortunately, that strategy did not yield good results in our experiments (data not shown). Another study was reported in FTA cards for the detection of several viral agents involved in the bovine respiratory complex, including bovine pestivirus, where viral RNA from respiratory tract swabs was stored for 14 days at temperatures between -27 and 46°C.

The article compared the RNA recovery from specimens in viral transport medium and FTA fixed samples, which proved 100% agreement¹³. Unlike the former, in this paper, the authors used a qPCR assay and carried out an elution protocol before RNA extraction. Finally, a third article describes the use of FTA cards to spot nasal swabs and blood samples from PI animals that were immediately processed by real time-PCR without storage. In this paper, the authors compared the mentioned sampling method to fresh samples, and observed good agreement between both, except for an increase in CT values in FTA samples compared with fresh samples, which was not considered significant given the high viral load in PI samples⁸.

Our study suggests that chromatography papers are a useful alternative for collecting and shipping blood samples for BDV diagnostic purposes, as a cost-effective method. However, for longer storage time (6 and 12 months), the recovery percentages decreased, although RNA was still detected, meaning that the biobanking purposes are not suitable for more than 12-month storage at any temperature with the protocols evaluated in this experiment. Nevertheless, the fact that pestiviral RNA was amplified in PI animal blood stored up to 12 months in chromatography paper, suggests that, although the recovery percentage was lower as the storage time progressed, it was still possible to recover RNA in blood samples from PI animals kept at refrigeration temperature.

FTA cards contain chelating agents and a free-radical trap designed to deal with atmospheric pollutants, thus protecting the entrapped nucleic acids for at least six years at room temperature¹. However, in our work, chromatography paper -which lacks these components-, yielded similar results. Moreover, for 12-month storage, drops spotted on chromatography paper showed slightly better results than on FTA cards, considering the number of positive repetitions.

When collecting samples, we recommend drying the drop spotted on filter papers at room temperature for at least 4 h (depending on ambient humidity), and sending them to the laboratory in individual zipped bags with a desiccant, such as silica gel. For shipping, refrigeration is not necessary.

The possibility of storing viral nucleic acids on filter papers constitutes a breakthrough in veterinary practice. The chance to collect samples directly on filter papers from herds where the disease is suspected could be very helpful, since in this kind of matrix, samples do not require immediate shipping or refrigeration and more importantly, they do not need to be processed instantly⁶. This methodology is mostly advantageous when performing molecular studies, since it guarantees nucleic acid stability at a wide range of temperatures and filter papers can be used with a variety of samples. Compared with the traditional method for collecting, shipping and storing samples, filter paper methodology represents a major improvement especially in large countries with poor infrastructure for shipping samples and also for those lacking laboratories which run this type of assays in every state or province.

An advantage of using chromatography paper as a method for collecting samples is its cost. Every sheet of paper can hold approximately 25 drops of sample, with a value of USD 66.55 per 100 sheets. In contrast, only four drops of sample can be spotted on every FTA card, with an approximate cost of USD 308 (depending on the provider) per 25 sheets.

Our study demonstrated that an economical filter paper, such as chromatography paper, is an appropriate alternative for collecting, shipping and storing blood samples for BVD diagnosis, being the results yielded as good as those using FTA specialized cards. The aforementioned is valid, as long as it is accompanied by a highly sensitive technique such as the RT-qPCR reported in this study.

Further studies using this methodology applied to other viruses are warranted, especially for diseases relevant to bovine health, which require mandatory control and are economically important for local production in order to improve sampling, shipping and laboratory processing conditions.

Funding

Financial support for the execution of this study was provided by Fundación ArgenINTA.

Conflict of interest

The authors declared no conflict of interest with respect to the research, authorship, and/or publication of this article.

Acknowledgements

The technical assistance of Maximiliano Jordan, Emiliano Nicodemo and Darío Malacari is greatly appreciated.

Appendix A. Supplementary data

Supplementary data associated with this article can be found, in the online version, at [doi:10.1016/j.ram.2022.10.002](https://doi.org/10.1016/j.ram.2022.10.002)

References

1. Ahmed HA, MacLeod ET, Hide G, Welburn SC, Picozzi K. The best practice for preparation of samples from FTA® cards for diagnosis of blood borne infections using African trypanosomes as a model system. *Parasit Vectors*. 2011;4:1–7.
2. Biswal JK, Subramaniam S, Ranjan R, Pattnaik B. Evaluation of FTA® card for the rescue of infectious foot-and-mouth disease virus by chemical transfection of extracted RNA in cultured cells. *Mol Cell Probes*. 2016;30:225–30.
3. Brock KV, Grooms DL, Ridpath J, Bolin SR. Changes in levels of viremia in cattle persistently infected with bovine viral diarrhoea virus. *J Vet Diagn Invest*. 1998;10:22–6.
4. Bustin SA, Benes V, Garson JA, Hellemans J, Huggett J, Kubista M, Mueller R, Nolan T, Pfaffl MW, Shipley GL. The MIQE guidelines: minimum information for publication of quantitative real-time PCR experiments. *Clin Chem*. 2009;55:1–12.
5. Cardona-Ospina JA, Villalba-Miranda MF, Palechor-Ocampo LA, Mancilla LI, Sepúlveda-Arias JC. A systematic review of FTA cards® as a tool for viral RNA preservation in fieldwork: are they safe and effective? *Prev Vet Med*. 2019;172:1–11.
6. Choi EH, Lee SK, Ihm C, Sohn YH. Rapid DNA extraction from dried blood spots on filter paper: potential applications in biobanking. *Osong Public Heal Res Perspect*. 2014;5:351–7.
7. Dubovi EJ. Laboratory diagnosis of bovine viral diarrhoea virus. *Biologicals*. 2013;41:8–13.

8. Foster C, Hoar B, Acvpm D, Hessman BE. Evaluation of FTA cards as a sample collection method for detection of BVDV in persistently infected beef cattle. *Bov Pract.* 2014;48:38–44.
9. Hall-Mendelin S, Hewitson GR, Genge D, Burtonclay PJ, De Jong AJ, Pyke AT, Van Den Hurk AF. FTA cards facilitate storage, shipment, and detection of arboviruses in infected aedes aegypti collected in adult mosquito traps. *Am J Trop Med Hyg.* 2017;96:1241–3.
10. Józwiak M, Wrosteck K, Domańska-Blicharz K, Olszewska-Tomczyk M, Śmietanka K, Minta Z. Application of FTA® Cards for detection and storage of avian influenza virus. *J Vet Res.* 2016;60:1–6.
11. Khodakaram-Tafti A, Farjanikish GH. Persistent bovine viral diarrhoea virus (BVDV) infection in cattle herds. *Iran J Vet Res.* 2017;18:154–63.
12. Lanyon SR, Hill FI, Reichel MP, Brownlie J. Bovine viral diarrhoea: pathogenesis and diagnosis. *Vet J.* 2014:201–9.
13. Liang X, Chigerwe M, Hietala SK, Crossley BM. Evaluation of Fast Technology Analysis (FTA) Cards as an improved method for specimen collection and shipment targeting viruses associated with Bovine Respiratory Disease Complex. *J Virol Methods.* 2014;202:69–72.
14. Linhares DCL, Rovira A, Torremorell M. Evaluation of Flinders Technology Associates cards for collection and transport of samples for detection of Porcine reproductive and respiratory syndrome virus by reverse transcription polymerase chain reaction. *J Vet Diagn Invest.* 2012;24:328–32.
15. MacLachlan NJ, Dubovi E. *Flaviviridae*. In: MacLachlan NJ, Dubovi E, editors. *Fenner's veterinary virology*. Academic Press; 2017. p. 525–45.
16. Meyling A, Mikél Jensen A. Transmission of bovine virus diarrhoea virus (BVDV) by artificial insemination (AI) with semen from a persistently-infected bull. *Vet Microbiol.* 1988;17:97–105.
17. Moennig V, Becher P. Control of bovine viral diarrhoea. *Pathogens.* 2018;7:1–12.
18. Monteiro FL, Cargnelutti JF, Martins B, Noll JG, Weiblen R, Flores EF. Detection of bovine pestiviruses in sera of beef calves by a RT-PCR based on a newly designed set of pan-bovine pestivirus primers. *J Vet Diagn Invest.* 2019;31:255–8.
19. Muñoz-Zanzi CA, Johnson WO, Thurmond MC, Hietala SK. Pooled-sample testing as a herd-screening tool for detection of bovine viral diarrhoea virus persistently infected cattle. *J Vet Diagn Invest.* 2000;12:195–203.
20. Newcomer BW. 75 years of bovine viral diarrhoea virus: current status and future applications of the use of directed antivirals. *Antiviral Res.* 2021;196:105205.
21. Nugroho W, Juniarti R, Silitonga P, Reichel MP, Irianingsih H, Wicaksono MS. The epidemiology and control of bovine viral diarrhoea virus in tropical Indonesian cattle. *Pathogens.* 2022;11:1–19.
22. Pfaffl MW. A new mathematical model for relative quantification in real-time RT-PCR. *Nucleic Acids Res.* 2001;29:2002–7.
23. Ririe KM, Rasmussen RP, Wittwer CT. Product differentiation by analysis of DNA melting curves during the polymerase chain reaction. *Anal Biochem.* 1997;245:154–60.
24. Robinson TL, Sutherland IA, Sutherland J. Validation of candidate bovine reference genes for use with real-time PCR. *Vet Immunol Immunopathol.* 2007;115:160–5.
25. Sakai T, Ishii A, Segawa T, Takagi Y, Kobayashi Y, Itou T. Establishing conditions for the storage and elution of rabies virus RNA using FTA® cards. *J Vet Med Sci.* 2015;77:461–5.
26. Sambrook J, Russell DW. Preparation and transformation of competent *E. coli* using calcium chloride. *Cold Spring Harb Protoc.* 2006;1:116–8.
27. Sarangi LN, Naveena T, Rana SK, Surendra KSNL, Reddy RVC, Bajjibabu P, Ponnanna NM, Sharma GK, Srinivasan VA. Evaluation of a specialized filter-paper matrix for transportation of extended bovine semen to screen for bovine herpesvirus-1 by real-time PCR. *J Virol Methods.* 2018;257:1–6.
28. Schirrmeier H, Strebelow G, Depner K, Hoffmann B, Beer M. Genetic and antigenic characterization of an atypical pestivirus isolate, a putative member of a novel pestivirus species. *J Gen Virol.* 2004;85:3647–52.
29. Smith DB, Meyers G, Bukh J, Gould EA, Monath T, Muerhoff AS, Pletnev A, Rico-Hesse R, Stapleton JT, Simmonds P, Becher P. Proposed revision to the taxonomy of the genus *Pestivirus*, family *Flaviviridae*. *J Gen Virol.* 2017;98:2106–12.
30. Smith RL, Sanderson MW, Walz PH, Givens MD. Sensitivity of polymerase chain reaction for detection of bovine viral diarrhoea virus in pooled serum samples and use of pooled polymerase chain reaction to determine prevalence of bovine viral diarrhoea virus in auction market cattle. *J Vet Diagn Invest.* 2008;20:75–8.
31. Spetter MJ, Louge Uriarte EL, Armendano JI, Álvarez I, Norero NS, Storani L, Pereyra SB, Verna AE, Odeón AC, González Altamiranda EA. Frequency of bovine viral diarrhoea virus (BVDV) in Argentinean bovine herds and comparison of diagnostic tests for BVDV detection in bovine serum samples: a preliminary study. *Brazilian J Microbiol.* 2021;52:467–75.
32. Spetter MJ, Louge Uriarte EL, Armendano JI, Morrell EL, Cantón GJ, Verna AE, Dorsch MA, Pereyra SB, Odeón AC, Saliki JT, González Altamiranda EA. Detection methods and characterization of bovine viral diarrhoea virus in aborted fetuses and neonatal calves over a 22-year period. *Braz J Microbiol.* 2020;51:2077–86.
33. Tam KI, Esona MD, Williams A, Ndze VN, Boula A, Bowen MD. Evaluation of BBL™ Sensi-Discs™ and FTA® cards as sampling devices for detection of rotavirus in stool samples. *J Virol Methods.* 2015;222:41–6.
34. Vilcek S, Strojny L, Durkovic B, Rossmanith W, Paton D. Storage of bovine viral diarrhoea virus samples on filter paper and detection of viral RNA by a RT-PCR method. *J Virol Methods.* 2001;92:19–22.
35. Weber MN, Mósena ACS, Simoes SVD, Almeida LL, Pessoa CRM, Budaszewski RF, Silva TR, Ridpath JF, Riet-Correa F, Driemeier D, Canal CW. Clinical presentation resembling mucosal disease associated with 'HoBi'-like pestivirus in a field outbreak. *Transbound Emerg Dis.* 2016;63:92–100.
36. Wernike K, Gethmann J, Schirrmeier H, Schroder R, Conraths FJ, Beer M. Six Years (2011–2016) of mandatory nationwide bovine viral diarrhoea control in Germany – a success story. *Pathogens.* 2017;6:1–8.